


# EPIDURAL MINI SET



## Tuohy Epidural Needle

- 16G or 18G needle with plastic stylet, fixed wing and 10mm interval depth markings
- Arrows marked on fixed wing assist needle placement

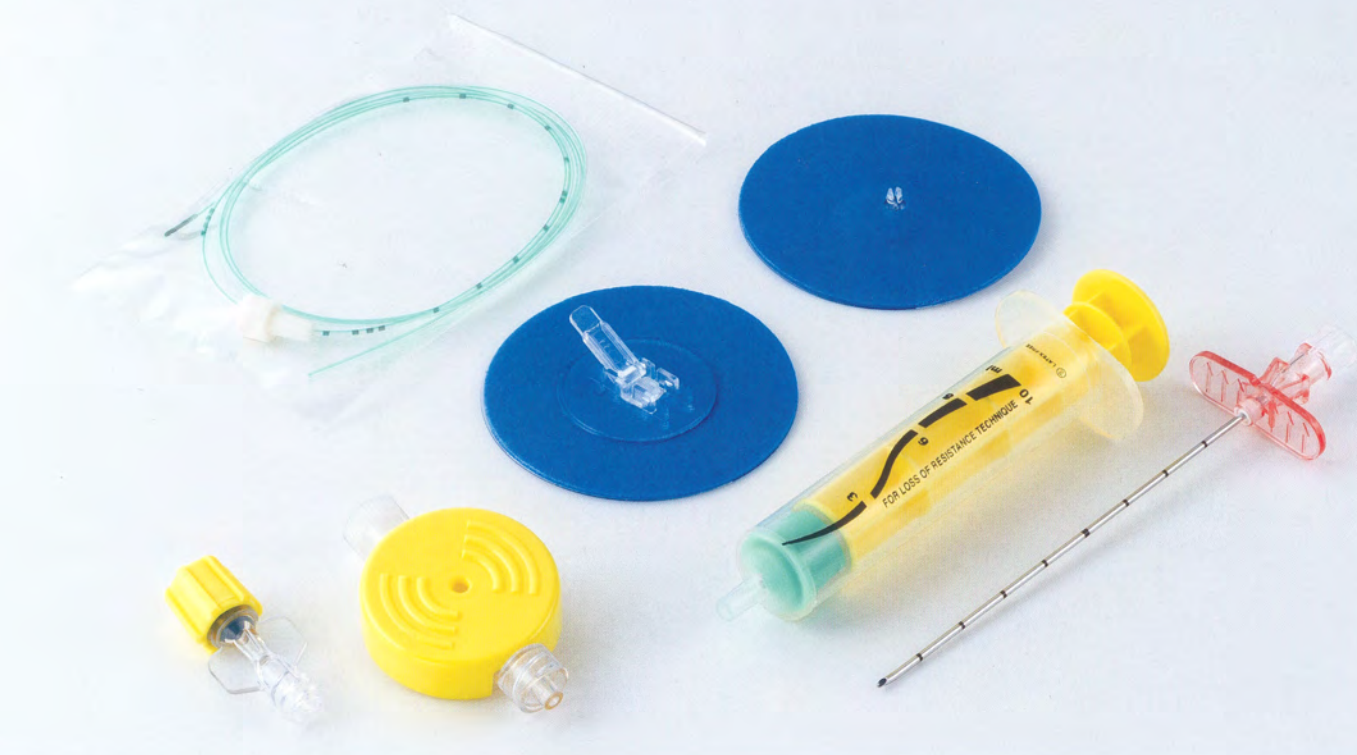

## Soft Tip Epidural Catheter

- Distance markings and radiopaque atraumatic closed soft tip prevent dura perforation
- 20G x 1000mm catheters with 3 lateral holes



## LOR Syringe (10ml)

- Plunger has an extremely low friction coefficient for ultra-smooth movement
- Suitable for use both with air and saline techniques




## Connector

- Ergonomically-designed Epidural Borst connector




## Filter

- 0.2 micron, 8 bar resistance epidural filter

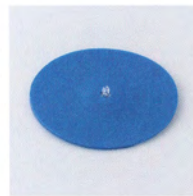

## Filter Fixing Device

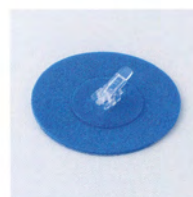

## Catheter Fixing Device

 **Dr. Japan Co., Ltd.**

1-1 Kagurazaka, Shinjuku-ku, Tokyo 162-0825 Japan  
<http://www.drjapan-jp.com>  
[info.overseas@drjapan-jp.com](mailto:info.overseas@drjapan-jp.com)

The Product specifications may change without notice.  
As of 4/2016 No. AN02-2